


VISION – Complete Suite of PC Software Options

Further enhance the productivity of your lab with the addition of powerful, easy-to-use VISION software to your UV1 spectrophotometer. Test setup, data collection, data manipulation, results presentation, regulatory compliance – all aspects of your analysis are easier with the VISION software series.

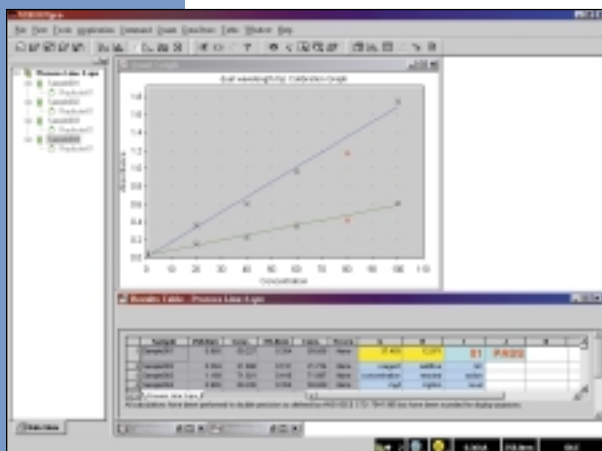

VISIONpro™

- A comprehensive UV-Vis software package with scanning, single wavelength, multiple wavelength, quantification and multi-component analysis supplied as standard functions
- Advanced data manipulation for easy method development and results calculation
- Advanced spreadsheet results calculation
- Spectral overlay and manipulation for easy comparison
- Easy-to-use and customizable interface
- Complete report flexibility including methods, graphs, data, Microsoft® Word files, Excel spreadsheets and graphics
- Multiple languages
- Requires Microsoft Windows® 98, 2000, NT or XP

VISIONsecurity™

- A complete software package with advanced security options and VISIONpro functionality
- Meets all security requirements of regulated and quality control laboratories
- Provides tools required to achieve 21 CFR Part 11 compliance
- Electronic Signatures for data and method traceability
- Audit trails for data and method manipulation history
- Requires either Microsoft Windows 2000/NT; or 98/2000/NT/XP with an NT or 2000 network

VISIONlife™

- Life science module for VISIONpro and VISIONsecurity
- Kinetics for intensity change with time studies
- Parallel rate measurement for advanced enzyme system studies
- Comprehensive manipulation and rate calculation options
- Enzyme kinetics and equation modeling

VISIONlite™

- Wavelength scanning with peak/valley identification, scan overlays, scale manipulation
- Rate measurement with enzyme activity calculation, data rescale/recalculation
- Fixed wavelength measurements at 1-31 wavelengths, absorbance ratio/difference with reference wavelength
- Quant measurement with tabular/graphical standard curve, linear & non-linear curve fitting, up to 20 standards
- Automatic data storage/recall, and data exporting to .csv file

Additional software packages available for the UV1 spectrophotometer include:

- Vision Auto™, used for automated analysis of multiple samples
- Vision Chroma™ for color measurement
- Enzymatic Food Analysis software
- Wine Analysis software

Highlights of VISION Software Functions

- *Single wavelength*
- *Multiple wavelength*
- *Quantification*
- *Multi-component analysis*
- *Scanning*
- *Kinetics*
- *Parallel rate*
- *Advanced enzyme modeling*
- *Report composer*
- *Simple password protection*
- *Multiple security levels*
- *Permission setting for software areas*
- *21 CFR Part 11 tools*
- *Electronic Signatures*
- *Color measurement*
- *Autosampler control*